



PROSonic Toothbrush

Where Style Meets Function™



PRO
SONIC

WHAT'S IN THE BOX

PROSonic Features:

- 1 PROSonic Brush with Rechargeable Lithium Ion Battery
- 1 Charger
- 1 Snap on Cap
- 2 Standard Brush Heads
- 1 USB cable
- 1 UL adapter

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE FIRST USE. STORE THIS MANUAL IN A SAFE PLACE FOR FUTURE REFERENCE!

Danger

To reduce the risk of electrocution:

1. Do not place or store product where it will be pulled or fall into a bath tub or sink.
2. Do not submerge the product in water or any other liquid.
3. Do not use while showering or bathing.
4. Do not reach for the product if it has fallen into water or other liquid. Unplug immediately if safe to do so.

Warnings

To reduce the risk of injury including burns, fire, or electrocution:

1. DO NOT plug into an outlet with a voltage other than specified on the bottom of the charger. This product is designed to operate within a range of 100 to 240 volts. Voltage converters DO NOT guarantee voltage compatibility.
2. DO NOT force the plug into an outlet if it does not easily fit. Discontinue use immediately.
3. DO NOT use the charger if it has a damaged cord.
4. DO NOT use the charger outdoors.
5. DO NOT use the charger if submerged in water.
6. DO NOT place any component of the product in the dishwasher for cleaning.
7. DO NOT use any attachments other than those recommended by the manufacturer.
8. Discontinue use if this product appears damaged in any way (brush head, handle, charger).
9. Keep cord away from all heat sources.
10. Do not use a brush head with damaged bristles. Crushed or bent bristles may break off while brushing. Replace the brush head every 3 months or sooner when signs of wear appear. Note that inappropriate cleaning may damage the bristles. (See Cleaning)

Warnings (Continued)

11. Use the PROSonic only for its intended use as described in this manual. This product was designed to only clean your teeth, gums and tongue. If pain or discomfort is experienced, discontinue use of this product and contact your dentist or physician as these may be signs of other health problems.
12. Persons (including children) with a lack of experience and knowledge or reduced mental, sensory or physical capabilities should not use this product, unless a person responsible for their safety has given proper instructions on the use of the product.
13. Children should be supervised to ensure proper usage and that they do not play with the toothbrush.
14. If your toothpaste includes peroxide, baking soda or other bicarbonate commonly found in whitening toothpastes, clean your brush head and handle thoroughly with water after each use. This will prevent plastic cracking due to prolonged exposure to these chemicals.

Medical Warnings

1. Consult your physician prior to using the PROSonic if you have any medical concerns.
2. If you have had oral or gum surgery in the previous two months, consult your dental professional prior to using this product.
3. If excessive bleeding occurs after using this product or bleeding continues to occur after one week of use contact your dental professional as this may be a sign of other health issues.
4. If you have a pacemaker or other implanted device contact your physician prior to use.
5. PROSonic has been tested and is compliant with safety standards for electromagnetic devices.

PRODUCT LIFE

Disposal

1. The rechargeable battery inside your PROSonic cannot be replaced, but should be removed from the handle for recycling.
2. If you have any questions about proper electronic waste disposal, contact your local waste management office.

WARRANTY

Guarantee & Service

POP SONIC guarantees the PRO4 Sonic for 1 year after the initial date of purchase. Defects due to faulty materials and/or workmanship will be repaired or replaced by POP SONIC when proof of purchase is provided during the warranty period. Contact our Customer Service Center at help@shoppopsonic.com

Warranty Exclusions

Items not covered under warranty:

1. Brush heads.
2. Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
3. Normal wear and tear, including scratches, abrasions, discoloration or fading.

Made by Pop Sonic, LLC
Minneapolis, MN 55344
www.shoppopsonic.com
Made in China

PRO SONIC

GETTING STARTED

Before using the POP SONIC PROSonic toothbrush for the first time, you will need to charge it. It will take approximately 4 hours to fully charge the PROSonic toothbrush for the first time or when the battery is completely depleted.

How to Charge

1. Plug the micro USB end of USB cord into the charge base, Plug the standard USB end of USB cord into EITHER the UL adapter OR another powered USB port. Plug the AC plug into a live outlet.
2. Place the PROSonic toothbrush on the charging stand as shown in the diagram. The battery indicator light on the toothbrush will turn on.
3. The battery indicator will automatically turn off when the charge cycle is complete.



Installing Brush Head

1. Insert the brush head onto the neck of the PROSonic toothbrush. Line up the guide mark on the back of the brush head with the start of the guide mark on the back of body of the brush as shown in the diagram
2. Rotate the brush head so that the guide marks on the back of the brush head and brush neck line up as shown in the diagram.

HOW TO USE

Brushing Instructions

1. Wet the brush head with water and add toothpaste.
2. Gently place the bristles of the brush on your teeth at an approximately 45° angle towards the gum line.

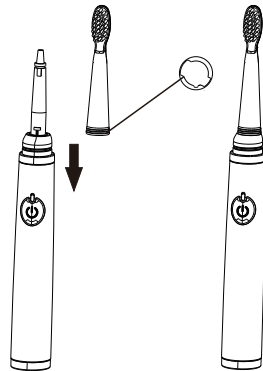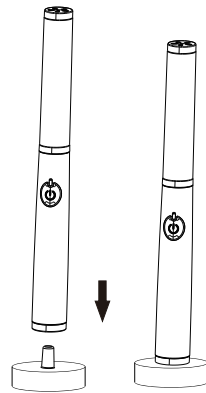

Dentists recommend that you replace PROSonic brush heads every 3 months to maintain optimal results.

Brushing Instructions (Continued)

3. Press the power button to activate the ultrasonic vibration.
4. Apply light pressure to let the PROSonic brush for you and achieve optimal results.
5. For the entire brushing cycle, slowly move the brush head gently across your teeth and gums in a back and forth motion. Be sure to allow the bristles to reach between your teeth and massage your gums.
6. Brush your teeth for 30 seconds on each of the quadrants of your mouth (4/30 method).
7. After brushing all four quadrants the sonic vibration will automatically stop after 2 minutes of continuous use, the dentist-recommended brushing time.
8. To stop at any time, press the power button. Note: The PROSonic is safe to use on braces and dental restorations (fillings, veneers, crowns, bridges, etc). When used on braces, the brush heads will wear out sooner.

Sonic Modes

On your first use, the PROSonic automatically starts with the default mode. Future use will use the previously selected mode. Select the appropriate mode by pressing the Mode Selection button. Paired with the proper included brush heads, the PROSonic toothbrush offers four modes for varying oral care needs.

Choose from the following:

- *Whiten (blue light indicator) : Removes surface stains*
- *Clean (green light indicator) : For exceptional everyday clean*
- *Sensitive (blue light flashing) : Gentle teeth and gum cleaning*
- *Massages (green light flashing) : Gently massages gums*

CLEANING & STORAGE

DO NOT CLEAN ANY PART OR COMPONENT IN THE DISHWASHER.

Brush Heads

1. Rinse the brush head and bristles after each use.
2. Shake to remove excess water.
3. Rinse the brush head connection weekly with warm water.

Toothbrush Cleaning

1. Remove the brush head and rinse the neck of the toothbrush with warm water.
2. Use a clean, damp cloth to wipe the entire surface of the toothbrush.

Storage

If you plan on not using the PROSonic for an extended period of time, unplug it, clean it, and store the entire unit in a cool, dry place and avoid direct sunlight.

TROUBLE SHOOTING

The PROSonic toothbrush won't turn on

You may need to recharge the PROSonic. Charge the PROSonic for 4 hours. (See Getting Started > How to Charge)

Tingling/tickling sensation when brushing

Due to the sonic vibrations, you may experience a tingling/tickling sensation when you first start using the PROSonic. This sensation will diminish over time as you become acclimated to the vibrations.

The brush head rattles against my teeth while brushing

Keep the brush head at a slight angle towards the gumline. Adjust the angle so the brush head doesn't touch the surface of your teeth.

The sonic vibrations seem weak

You may need to fully charge the PROSonic. You may need to replace your brush head. Also note that each mode has a distinct vibration pattern