

Product Details

Product Name: Prosonic 2.0
Input: 5V = 1.5W
Product material: ABS/TPE/Dupont Nylon
Bristle diameter: 0.127mm, 0.152mm soft
Toothbrush size: 26.9 X 26.4 x 247mm
Charging Time: ~4 hours
Working Temperature: 0-50°C
Origin: China

Warranty

Pop Sonic warrants the Pro 2.0 for a period of 1 year after the original date of purchase against defects due to faulty workmanship or materials. For warranty claims, you must register your device and log in to your account at www.shoppopsonic.com, then select the option to make a warranty claim. Our customer service department will contact you to initiate the warranty procedure.

This warranty covers working parts. It does not cover damage due to normal usage or accidents such as dropping the device.

Other exclusions include misuse or neglect to the device. This warranty is void if any other charging cables besides the one you received is used with your device.

For any warranty claims, proof of purchase must be provided. For this reason, please keep your receipt for the duration of the warranty period.

Shipping costs are non-refundable.

For questions about your product, please call Pop Sonic customer services at (800)772-9174 or email help@shoppopsonic.com

Made by Pop Sonic
Minneapolis, MN 55344 USA

MADE IN CHINA

IMPORTANT SAFETY NOTES

1. Users new to electric toothbrushes may experience slight numbing or other sensations after the first few uses.
2. Do not disassemble, repair, or alter the product without authorization
3. Do not break or drop the product. IF the product is defective, please discontinue use.
4. Please verify that your local voltage is compliant with the marked voltage before connecting to power.
5. This product is compatible with ProSonic brush heads only. Replacement heads can be purchased at www.shoppopsonic.com
6. If you’ve had oral or gum surgery in the past two months, please consult a dentist before using Pro 2.0.
7. If you experience bleeding after repeated use and/or the bleeding continues for more than seven days, please consult your dentist.
8. This product is designed for cleaning your teeth, gums, and tongue, only.
9. Children should be supervised to ensure proper usage. This product is not a toy.
10. Avoid placing in high temperature conditions or in direct sunlight.



Troubleshooting

If you’re experiencing weaker-than-normal vibrations or your brush time is cut short, try applying less pressure while you brush or verify that the brush is charging correctly.


If the product is not operating at all, contact customer service.



PRO 2.0

The ultimate clean just got better.
Sonic Toothbrush & Brush heads

PRO 2.0

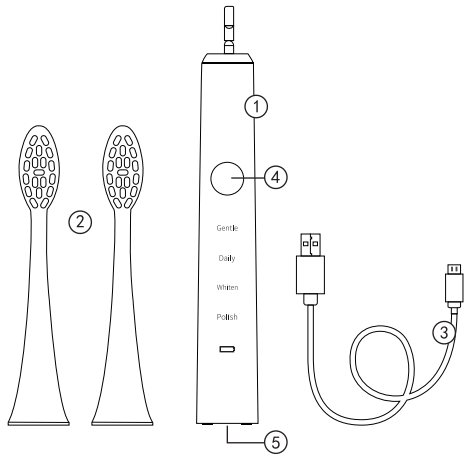

What’s Included

1. Toothbrush handle
2. Two brush heads
3. UBS charging cable

Components

4. Multifunction button & indicator light
5. Charging port


Features:

- 40000 Vibrations per minute
- 4 Cleaning modes (noted by indicator light): gently, daily, whiten, polish
- Anti-skid ergonomic design
- IPX7 waterproof rating
- Safe & durable lithium battery runs for up to 21 days on a full charge
- Low battery alerts
- USB charger
- Smart mode memory
- Two-minute auto stop function
- 30 Second vibrating timer to indicate how long to spend on each section of the mouth

Assembly & Disassembly

Gently push the brush head onto the handle until it locks into place.

Simply pull it off when it's time to replace your brush head. Replacement brush heads can be found at www.shoppopsonic.com

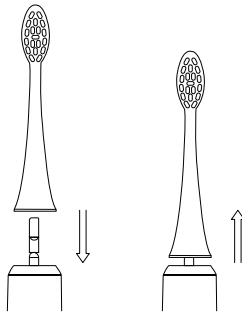

Using your Toothbrush

1. Using an AC outlet, plug the provided cable into the brush to begin charging.
2. Attach the brush head to the handle, and lock in place.
3. Wet the bristles and apply toothpaste.
4. Gently hold the brush head to your teeth.
5. Press the multifunction button and choose mode according to personal preference.
6. Gently move the brush head across each section of the teeth for 10 seconds, including: front, back, and chewing surfaces
7. Move to the next section when the auto reminder pauses briefly every 30 seconds
8. After a two-minute period, the brush will turn off.

Charging your Toothbrush

1. When battery capacity is low, the charging indicator light will flash 10 times.
2. Connect the toothbrush to the USB cable, then connect the plug to a power source to begin charging. The red charging indicator light will illuminate
3. When the battery is fully charge, the indicator light and charging unit will shut off automatically to avoid overcharging.

- At first use or when the battery has been dead for a period of time, it's common for the indicator light not to illuminate until it's been charging for five minutes or more.
- 6.

Lithium Battery Replacement and Discard

Do not attempt to remove lithium ion battery. These types of batteries carry inherent dangers if mishandled and can explode, catch fire, or cause burns if disassembled or are damaged.

Lithium batteries are subject to local disposal regulations and should not be discarded in regular trash. Check your local laws and regulations for disposal information.

Lithium batteries can be recycled. For recycling information, contact your local recycling service or visit online recycling resources such as www.batteryrecycling.com or www.rbrc.com.

