




USBSonic Toothbrush

Where Style Meets Function

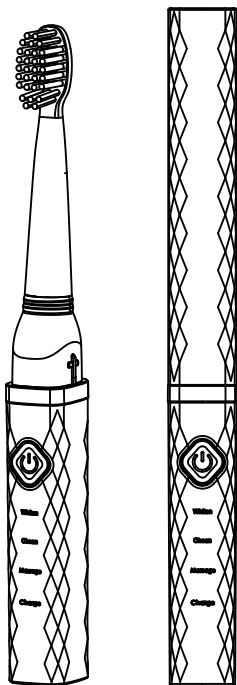

USB
SONIC

WHAT'S IN THE BOX

USBSonic Features:

- 1 USBSonic Brush with Rechargeable Lithium Ion Battery
- 1 USB Cord
- 1 Snap on Cap
- 2 Standard Brush Heads

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE FIRST USE. STORE THIS MANUAL IN A SAFE PLACE FOR FUTURE REFERENCE!

Danger

To reduce the risk of electrocution:

1. Do not place or store product where it will be pulled or fall into a bath tub or sink.
2. Do not submerge the product in water or any other liquid.
3. Do not use while showering or bathing.
4. Do not reach for the product if it has fallen into water or other liquid. Unplug immediately if safe to do so.

Warnings

To reduce the risk of injury including burns, fire, or electrocution:

1. DO NOT place any component of the product in the dishwasher for cleaning.
2. DO NOT use any attachments other than those recommended by the manufacturer.
3. Discontinue use if this product appears damaged in any way (brush head, handle).
4. Keep cord away from all heat sources.
5. Do not use a brush head with damaged bristles. Crushed or bent bristles may break off while brushing. Replace the brush head every 3 months or sooner when signs of wear appear. Note that inappropriate cleaning may damage the bristles. (See Cleaning)

Warnings (Continued)

6. Use the USBSonic only for its intended use as described in this manual. This product was designed to only clean your teeth, gums and tongue. If pain or discomfort is experienced, discontinue use of this product and contact your dentist or physician as these may be signs of other health problems.
7. Persons (including children) with a lack of experience and knowledge or reduced mental, sensory or physical capabilities should not use this product, unless a person responsible for their safety has given proper instructions on the use of the product.
8. Children should be supervised to ensure proper usage and that they do not play with the toothbrush.
9. If your toothpaste includes peroxide, baking soda or other bicarbonate commonly found in whitening toothpastes, clean your brush head and handle thoroughly with water after each use. This will prevent plastic cracking due to prolonged exposure to these chemicals.

Medical Warnings

1. Consult your physician prior to using the USBSonic if you have any medical concerns.
2. If you have had oral or gum surgery in the previous two months, consult your dental professional prior to using this product.
3. If excessive bleeding occurs after using this product or bleeding continues to occur after one week of use contact your dental professional as this may be a sign of other health issues.
4. If you have a pacemaker or other implanted device contact your physician prior to use.
5. USBSonic has been tested and is compliant with safety standards for electromagnetic devices.

PRODUCT LIFE

Disposal

1. The rechargeable battery inside your USBSonic cannot be replaced, but should be removed from the handle for recycling.
2. If you have any questions about proper electronic waste disposal, contact your local waste management office.

WARRANTY

Guarantee & Service

POP SONIC guarantees the USBSonic for 1 year after the initial date of purchase. Defects due to faulty materials and/or workmanship will be repaired or replaced by POP SONIC when proof of purchase is provided during the warranty period. Contact our Customer Service Center at www.shoppopsonic.com

Warranty Exclusions

Items not covered under warranty:

1. Brush heads.
2. Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
3. Normal wear and tear, including scratches, abrasions, discoloration or fading.

USB SONIC

GETTING STARTED

Before using the POP SONIC USB Sonic toothbrush for the first time, you will need to charge it. It will take approximately 8 hours to fully charge the USB Sonic toothbrush for the first time or when the battery is completely depleted.



How to Charge

1. Open the USB cover at the bottom of the toothbrush
2. Plug the micro USB end of the USB cord into the end of the toothbrush
3. Plug the standard USB end of the USB cord into EITHER the provided USB adapter OR another powered USB port
4. Plug the adapter into a live outlet.

Installing Brush Head

1. Insert the brush head onto the neck of the USB Sonic toothbrush. Line up the guide mark on the back of the brush head with the start of the guide mark on the back of body of the brush as shown in the diagram
2. Rotate the brush head so that the guide marks on the back of the brush head and brush neck line up as shown in the diagram.

Dentists recommend that you replace USB Sonic brush heads every 3 months to maintain optimal results.

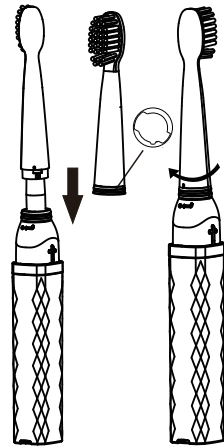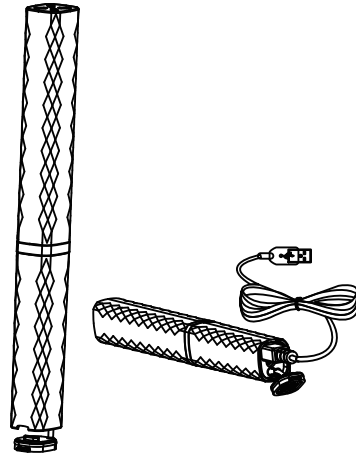

HOW TO USE Brushing Instructions

1. Wet the brush head with water and add toothpaste.
2. Gently place the bristles of the brush on your teeth at an approximately 45° angle towards the gum line.
3. Press the power button to activate the sonic vibration.

Brushing Instructions (Continued)

4. Apply light pressure to let the USB Sonic brush for you and achieve optimal results.
5. For the entire brushing cycle, slowly move the brush head gently across your teeth and gums in a back and forth motion. Be sure to allow the bristles to reach between your teeth and massage your gums.
6. Brush your teeth for 30 seconds on each of the quadrants of your mouth (4/30 method).
7. After brushing all four quadrants the sonic vibration will automatically stop after 2 minutes of continuous use, the dentist-recommended brushing time.
8. To stop at any time, press the power button. Note: The USB Sonic is safe to use on braces and dental restorations (fillings, veneers, crowns, bridges, etc). When used on braces, the brush heads will wear out sooner.

Sonic Modes

On your first use, the USB Sonic automatically starts with the default mode. Future use will use the previously selected mode. Select the appropriate mode by pressing the Mode Selection button. Paired with the proper included brush heads, the USB Sonic toothbrush offers three modes for varying oral care needs.

Choose from the following:

- Whiten
Removes surface stains
- Clean
For exceptional everyday clean
- Massage
Gently massages gums

CLEANING & STORAGE

DO NOT CLEAN ANY PART OR COMPONENT IN THE DISHWASHER.

Brush Heads

1. Rinse the brush head and bristles after each use.
2. Shake to remove excess water.
3. Rinse the brush head connection weekly with warm water.

Toothbrush Cleaning

1. Remove the brush head and rinse the neck of the toothbrush with warm water.
2. Use a clean, damp cloth to wipe the entire surface of the toothbrush.

Storage

If you plan on not using the USB Sonic for an extended period of time, unplug it, clean it, and store the entire unit in a cool, dry place and avoid direct sunlight.

TROUBLE SHOOTING

The USB Sonic toothbrush won't turn on

You may need to recharge the USB Sonic. Charge the USB Sonic for 8 hours. (See Getting Started > How to Charge)

Tingling/tickling sensation when brushing

Due to the sonic vibrations, you may experience a tingling/tickling sensation when you first start using the USB Sonic. This sensation will diminish over time as you become acclimated to the vibrations.

The brush head rattles against my teeth while brushing

Keep the brush head at a slight angle towards the gumline. Adjust the angle so the brush head doesn't touch the surface of your teeth.

The sonic vibrations seem weak

You may need to fully charge the USB Sonic. You may need to replace your brush head. Also note that each mode has a distinct vibration pattern